



© Primal Pictures Ltd, 2000



This model contains these conductive anatomical points:


- * Trigger finger
- * Carpal tunnel space
- * First carpal metacarpal joint
- * DeQuervain's disease
- * Radioulnar & carpal joint

Locate the appropriate spot for injection (i.e. joint space bursa).

Insert needle into location. The audio feedback box will emit a tone when the tip of the needle has made contact with the correct area.

Use the dry erase marker to identify injection targets.

Hand/Wrist Conductive Injection Model Item # 1511-5

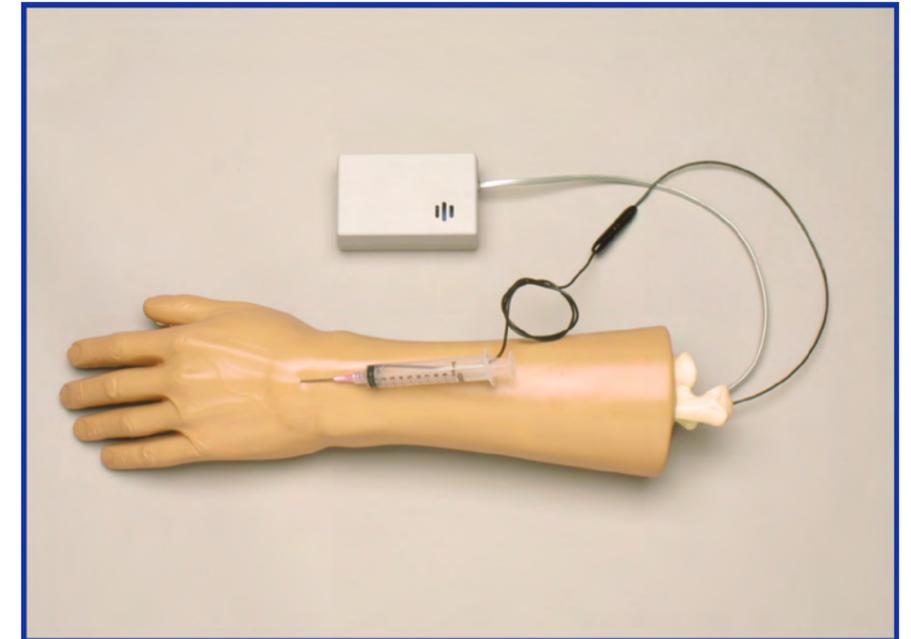

Thank you for purchasing the Sawbones conductive injection model. This injection trainer can be used to teach or demonstrate needle placement for aspiration and joint injection. With the correct needle insertion, an audible tone will sound indicating proper placement of the chosen area. Palpation of the major anatomical landmarks may also be performed. Additional conductive structures may be added upon request. With a little care and maintenance your new model should provide you with many years of service.

This instructional brochure will take you through:

- * Components
- * Set up instructions
- * Operation and anatomical features
- * Care and maintenance
- * Available replacement parts

SAWBONES®

Should you have any questions, please contact:


SAWBONES
10221 SW 188th Street
Vashon, WA 98070
(206) 463-5551
info@sawbones.com

SAWBONES EUROPE AB
Krossverksgatan 3
216 16 Malmö
Sweden
+46 40 650 70 00
info@sawbones.se

www.sawbones.com

SAWBONES®


COMPONENTS

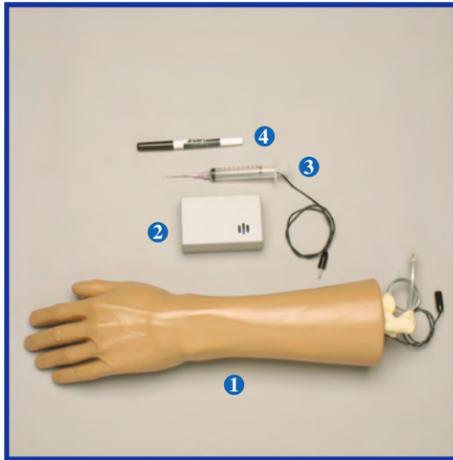

The following components are included:

- 1 Soft tissue injection hand/wrist
- 2 Audio feedback box with leads
- 3 Needle assembly with lead
- 4 Dry erase marker

SET UP

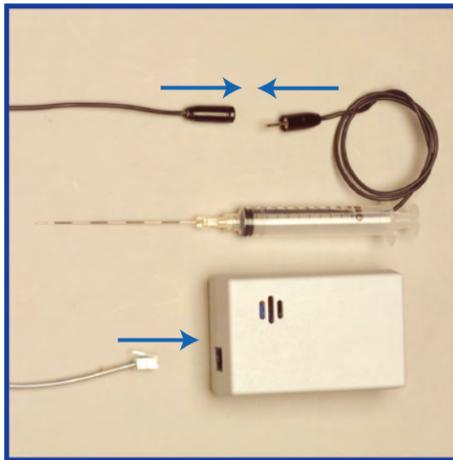

- 1 Lay out all parts and inspect for shipping damage.
- 2 Insert wire lead from syringe into the black wire.
- 3 Insert male phone connector into audio feedback box.



CARE & MAINTENANCE




IMPORTANT


To maintain audio feedback, do not remove plunger assembly or needle from the syringe. This model is NOT intended for injection of fluids.

Your model can be wiped down with any mild cleaner. Use only the dry erase marker on the model. Standard ink can not be removed. Store in a dry, dark place out of direct contact with sunlight.

REPLACEMENT PARTS

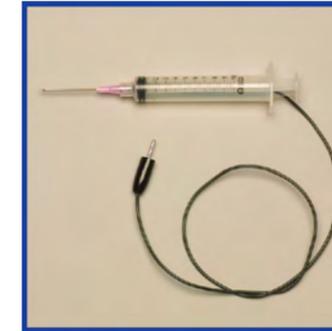

Needle assembly with lead
Item # 1620-1

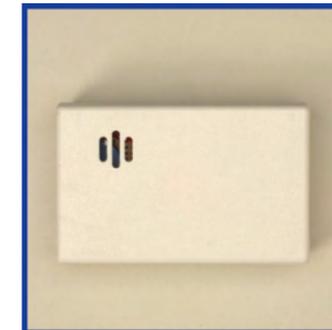

Conductive audio feedback box
Item # 1620



Dry erase marker
Item # 1903